




Silicone Lite Foam Border & Non Border

3

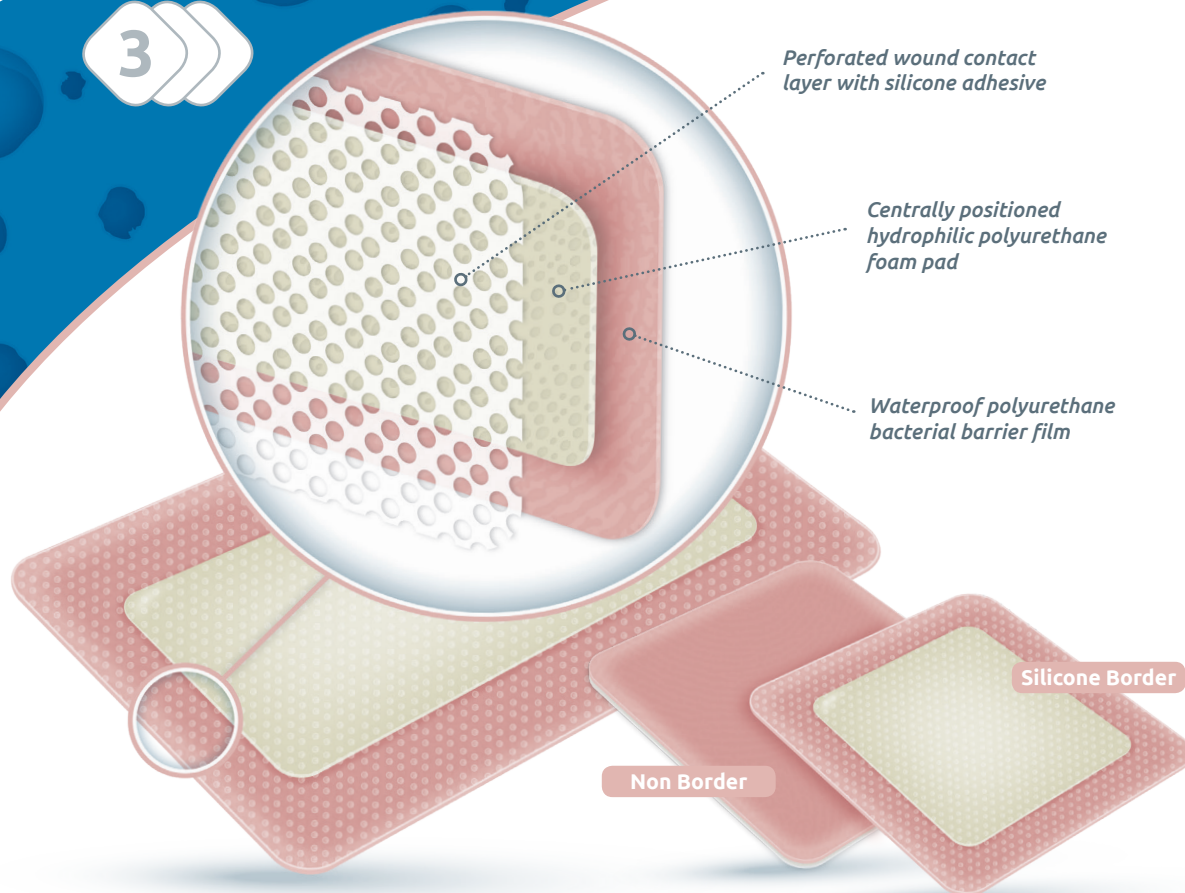

Perforated wound contact layer with silicone adhesive

Centrally positioned hydrophilic polyurethane foam pad

Waterproof polyurethane bacterial barrier film

Non Border

Silicone Border



Advanced Medical Solutions Ltd

Thin, soft and discreet dressing which manages non to low exuding wounds

Key Features

- ✓ **Soft silicone atraumatic adhesive**
Dressing can be repositioned during application or lifted during wear time for observations. Silicone wound contact layer aids atraumatic removal.
- ✓ **Thin, soft and discreet hydrophilic foam pad**
Dressing is able to conform to different anatomical contours.
- ✓ **Waterproof and bacterial barrier**
Protects both the patients and healthcare professional from cross contamination. Keeps bacteria and other contaminants away from the wound. The low friction backing enables the dressing to be used under compression bandaging.

Artist's impression

The hydrophilic foam has a flexible and breathable backing and a film which ensures that moisture can be transpired through the top of the dressing and works with the unique ability of the foam to absorb exudate.

The wound contact side of the dressing has a clear silicone perforated film which aids atraumatic removal, and prevents ingress of granulation tissue into the dressing. The perforations in the wound contact layer allows the uptake of exudate to minimise maceration and damage to the surrounding periwound skin.



Key benefits

- ✓ **Soft and conformable foam pad**
Easy application. Conformable to anatomical contours.
- ✓ **Minimal epidermal stripping or pain on removal**
Dressing can be lifted and repositioned during wear time.


Product selection



Product Selection			
📐 Border		📐 Non Border	
5cm x 5cm	2" x 2"	7.5cm x 7.5cm	3" x 3"
7.5cm x 7.5cm	3" x 3"	10cm x 10cm	4" x 4"
10cm x 10cm	4" x 4"	15cm x 15cm	6" x 6"
15cm x 15cm	6" x 6"	10cm x 20cm	4" x 8"
📐 Operative Border			
5.5cm x 12cm	2" x 4.7"		
10cm x 15cm	4" x 6"		
10cm x 20cm	4" x 8"		
10cm x 25cm	4" x 10"		
10cm x 30cm	4" x 12"		

Performance

Total Fluid Handling¹

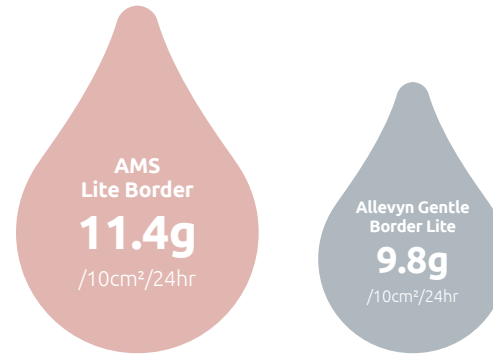

Indications

For use in the management of non to low exuding, chronic and acute wounds throughout the healing process. The dressing is indicated for use on the following wounds:

- ▶ Pressure ulcers
- ▶ Diabetic ulcers
- ▶ Leg & foot ulcer
- ▶ Post-operative surgical wounds
- ▶ Trauma wounds (including lacerations, abrasions, skin tears and blisters)
- ▶ Superficial and partial thickness burns
- ▶ Radiation damaged skin

Contact

Advanced Medical Solutions Ltd
Premier Park
33 Road One
Winsford Industrial Estate
Winsford
Cheshire
CW7 3RT
UK

Tel: +44(0) 1606 863 500
Email: sales@admedsol.com

www.admedsol.com



Helping you help others

Advanced Medical Solutions excels in the research and development of foam dressings, through continuous product development and by working closely within the market place we meet the needs of the clinician, patient and our commercial partners.

References: 1. AMS data on file – LD033-17 and LD068-17
Allevyn Gentle Border Lite is a registered trademark of Smith and Nephew.