



Enhanced performance from a next generation silver alginate dressing

Silver IV Performance

- ▶ High gelling properties and absorbency have not compromised the wet strength
- Higher absorbency can result in reduced dressing change frequency, minimising clinician time and disruption to the patient¹

Key Features

High silver content¹

Improved antimicrobial performance.

Fast initial silver elution²

Quickly establishes an effective level of antimicrobial activity.

Sustained silver release²

On-going silver ion release allowing antimicrobial effects for up to 7 days.

High wet strength¹

Can be removed intact minimising discomfort at dressing change.

Artist's impression

Silver IV Antimicrobial wound dressing is a sterile, non-woven pad consisting of a high M (mannuronic acid) calcium alginate and carboxymethylcellulose (CMC).

A coloured ionic silver salt releases silver ions in the presence of wound fluid. As fluid is absorbed, the alginate forms a soft gel which assists in maintaining a moist environment for optimum wound healing and allows intact removal. Silver ions released in the presence of wound exudate are an effective antimicrobial agent for up to 7 days, based on *in-vitro* testing, against a broad spectrum of microorganisms frequently associated with bacterial colonisation and infection of wounds.





AMS Silver IV Antimicrobial wound dressing is effective against the following microorganisms, based on in-vitro³ testing:

- Staphylococcus epidermidis including MRSA
- Streptococcus pyogenes including MRSE
- Enterococcus faecalis (VRE)
- Staphylococcus aureus
- Pseudomonas aeruginosa
- Escherichia coli
- Candida albicans

Key benefits

Gelling fibres

As wound fluid is absorbed, the fibres form a soft gel which follows the contours of the wound bed.

Reduces odour

Due to the antimicrobial properties the dressing can help reduce the odour of infected wounds¹.

Indications²


Silver IV Antimicrobial wound dressing is indicated for the following moderately to heavily exuding, partial to full thickness wounds:

- Post operative wounds
- Trauma wounds (dermal lesions, trauma injuries or incisions)
- Leg ulcers

- Pressure ulcers
- Diabetic ulcers
- ▶ Graft and donor sites
- Cavity wounds
- Superficial and partial thickness burns

24% faster Silver Elution (PPM) than Aquacel Ag Extra

Within the first 12 hours¹

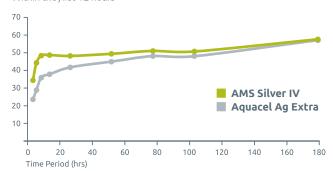

2.7cm x 3

Product selection

Product Selection

5cm x 5cm

10cm x 10cm

10cm x 12cm

10cm x 20cm

15cm x 15cm

20cm x 30cm



		_			
32cm					
	 • • •		 	•	

Contact

Advanced Medical Solutions Ltd Premier Park 33 Road One Winsford Industrial Estate Winsford Cheshire CW7 3RT UK

Tel: +44(0) 1606 863 500 Email: sales@admedsol.com

www.admedsol.com



Helping you help others

Advanced Medical Solutions is a world leader in the research, development, manufacture and supply of alginates, co-spun alginates and other proprietary products to global advanced wound care partners.